# EMOTECHE *BLOOD BANKS*

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## **IMPORTANT WARNING**

The high level of technology of the products made by the manufacturer occasionally requires components to be adapted from parts normally found on the market in their standard commercial version. Therefore, if it should be ascertained that any parts have been replaced with other components that are not original, i.e. that have not been supplied and guaranteed directly by the MANUFACTURER, or that any unauthorised technical interventions have been carried out, the guarantee, if any is in force, shall be considered no longer valid and the MANUFACTURER himself shall no longer be held responsible towards the purchaser and towards third parties.

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## SUMMARY

<b>1</b>	<b>SUMMARY OF RATING PLATE DATA</b>	<b>2</b>
<b>2</b>	<b>WARNINGS</b>	<b>3</b>
2.1	GENERAL WARNINGS	3
2.2	WARNINGS FOR TRANSPORT AND HANDLING	4
2.3	WARNINGS FOR INSTALLATION	4
2.4	WARNINGS FOR PERSONNEL IN CHARGE OF THE MACHINE	4
2.5	WARNINGS FOR MAINTENANCE	5
2.6	EXPLANATION OF SYMBOLS	5
<b>3</b>	<b>TECHNICAL SPECIFICATIONS</b>	<b>7</b>
3.1	TECHNICAL DATA	7
3.2	FUSES	8
3.3	ENVIRONMENTAL CONDITIONS	8
3.4	TECHNICAL LAY-OUT DIAGRAMS	9
<b>4</b>	<b>HANDLING AND REMOVAL OF PACKAGING</b>	<b>10</b>
4.1	PERSONNEL REQUISITES	10
4.2	MACHINE CONDITIONS	10
4.3	EQUIPMENT NEEDED FOR HANDLING	10
4.4	HOW TO RAISE AND MOVE WITH A BRIDGE CRANE	11
4.5	REMOVAL OF PACKAGING	11
4.5.1	HOW TO REMOVE THE OUTER STRAPS	12
4.5.2	HOW TO REMOVE THE CARDBOARD CASING	12
4.5.3	REMOVAL OF PLASTIC MATERIALS	12
4.5.4	HOW TO REMOVE THE WOODEN SIDE BLOCKS	12
4.5.5	HOW TO REMOVE THE PALLET	13
4.6	HOW TO BLOCK THE MACHINE	13
<b>5</b>	<b>DESCRIPTION OF THE SYSTEM</b>	<b>14</b>
5.1	GENERAL VIEW	14
5.2	INTERNAL LIGHT	15
5.3	CONTROL SYSTEM	15
<b>6</b>	<b>INSTALLATION</b>	<b>16</b>
6.1	POSITIONING OF THE MACHINE	16
6.2	ELECTRICAL WIRING	17
6.2.1	HOW TO CONNECT TO THE ELECTRICAL MAINS SUPPLY	17
6.2.2	SUPPLY CABLE	18
6.2.3	INSERTION OF REMOTE ALARM	18
6.2.4	SERIAL INTERFACE	18
<b>7</b>	<b>USE FORESEEN BY THE MANUFACTURER</b>	<b>19</b>
7.1	AIM OF THE MACHINES	19
7.2	INTENDED USE OF THE MACHINES	19
7.3	OPERATOR	19
7.4	USE LIMITS	19
7.5	USE OF PROTECTIVE CLOTHING	19
<b>8</b>	<b>START-UP</b>	<b>20</b>
8.1	E2003 SYSTEM	20
<b>9</b>	<b>USE</b>	<b>20</b>
9.1	PRODUCT LOADING	21
9.2	TEMPERATURE CHART RECORDER	22
9.2.1	USE OF THE RECORDER	22
9.3	HOW TO INSERT AND CHANGE THE DRAWERS (OPTIONAL)	24
9.4	CONNECTION OF OPTIONAL PT100 PROBE (OPTIONAL)	24
9.5	DATA DOWNLOADING VIA SERIAL PORT (OPTIONAL)	25
<b>10</b>	<b>SAFETY DEVICES- CHECK AND SET-UP</b>	<b>25</b>
<b>11</b>	<b>MAINTENANCE</b>	<b>26</b>
11.1	ALARM SYSTEM CONTROL	26
11.2	DEFROSTING	26
11.3	HOW TO CLEAN THE STRUCTURE	26
11.3.1	HOW TO CLEAN THE GASKETS	27
11.3.2	HOW TO CLEAN THE CONDENSER	27
11.4	COOLING SYSTEM MAINTENANCE	27
11.5	MAINTENANCE OF THE ELECTRICAL-ELECTRONIC PART	27
<b>12</b>	<b>TROUBLESHOOTING</b>	<b>28</b>
12.1	HOW TO REPLACE THE LAMP	29
<b>13</b>	<b>REMOVAL FROM INSTALLATION SITE</b>	<b>30</b>
13.1	DISASSEMBLY	30
13.2	SCRAPPING	30
<b>14</b>	<b>DIAGRAMS</b>	<b>31</b>
14.1	REFRIGERATOR DRAWING	31
14.2	ELECTRIC DIAGRAM	32

  	
laboratory and biomedical division	
MORE THAN TECHNOLOGY	
Model	<input type="text"/>
Serial No	<input type="text"/>
Year	<input type="text"/>
Weight kg	<input type="text"/>
°C	<input type="text"/>
RH%	<input type="text"/>
IP	<input type="text"/>
POWER SUPPLY	V~ <input type="text"/> Hz <input type="text"/> P Kw <input type="text"/> A <input type="text"/> Class <input type="text"/>
REFRIGERANT PROPERTIES	n. 1 <input type="text"/> Pmax kPa <input type="text"/> n. 2 <input type="text"/> Pmax kPa <input type="text"/>
COUNTRY OF ORIGIN: ITALY	

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## GB 1 SUMMARY OF RATING PLATE DATA

- Find the special rating plate on the machine in order to check the technical data.
- Check the machine model and supply voltage before carrying out any operation whatsoever.
- If you find any discrepancies, contact the manufacturer or your supplier immediately.

1	EC compliancy mark	9	Electrical protection grade
2	Manufacturer's address	10	Supply voltage
3	Model	11	Absorbed power
4	Series N°	12	Electrical absorption
5	Year of manufacture	13	Frequency
6	Temperature	14	Isolation class
7	Relative humidity	15-16	Type of refrigerant
8	Weight	17-18	Pressure of refrigerant

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### 2.1 GENERAL WARNINGS

- Do not carry out any operations or manoeuvres unless you are absolutely certain of their effect; if in doubt, contact your nearest technical assistance service or the manufacturer himself.
- The manufacturer will not be held responsible for damage to the machine or to objects in the following cases:
  - improper use of the machine
  - use of unsuitable personnel
  - incorrect assembly and installation
  - defects in the plant systems
  - unauthorized modifications or operations to the machine
  - use of spare parts that are not original pieces
  - failure to comply with the norms given in this handbook
  - exceptional events
- This instruction handbook has been designed for the following personnel:
  - Personnel in charge of transport, handling and removal of packaging
  - Personnel in charge of the preparation of the plant systems and installation site
  - Installers
  - Personnel in charge of using the machine
  - Personnel in charge of maintenance
- The instruction handbook indicates the use foreseen by the manufacturer and cannot ever replace adequate experience of the operator. It can only be used as a reminder of the main operations to be carried out.
- The instruction handbook should be kept carefully and should also be within easy reach for reference. If necessary, photocopy the pages concerned directly with machine use.

The handbook should last at least the life-time of the machine itself.

- The instruction handbook gives technical information on how the machine is manufactured at the present time; the manufacturer reserves the right to carry out any modifications he deems necessary to the machines and to the instruction handbooks, without prior notice or replacement.

If you lose or destroy the handbook, you may ask for a copy. Please give the rating plate data (see chap. 1) in your request.

**This machine is not referred to in enclosure IV of the community norm on machines, and therefore the procedure in article 8, paragraph 2, letter A in the Community norms for machines 89/392/CEE and 91/368/EEC has been applied. Article 8, paragraph 2, letter A of the Community norms for machines 89/392/CEE and 91/368/EEC obliges the manufacturer to provide the file foreseen by the enclosure of the afore-mentioned laws, and to store it in the technical archives of the company Angelantoni Industrie spa, località Cimacolle, Massa Martana, Perugia.**

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### 2.2 WARNINGS FOR TRANSPORT AND HANDLING


- This symbol, placed on each packaging, indicates the weight of each package. Always check that the tools and machines to handle and transport the machine are adequate.
- Always keep the machine in an upright position. If the machine should accidentally turn upside-down or on its side, do not switch it on. Put it in the correct position and contact the manufacturer.

### 2.3 WARNINGS FOR INSTALLATION

- Installation should always be carried out by specialized personnel.
- Carefully follow the instructions on how to prepare the plant systems before installing the machine.
- When the installation site is being prepared, bear in mind the space and work conditions of the personnel in charge of the machine so as to reduce to a minimum noise, fatigue, discomfort and anything else which may have a negative influence on the staff.
- When designing the installation site, remember to leave sufficient space for control, maintenance, cleaning, and removal of production waste material.
- Make sure that the work site is adequately lit so that personnel can work in optimum conditions.
- When designing the installation site, please refer to the norms in force and in particular:
  - set up all the firefighting and safety devices.

### 2.4 WARNINGS FOR PERSONNEL IN CHARGE OF THE MACHINE

- The machine may only be used by personnel who have read the rules described in this handbook.
- To open the chamber when its temperature is very different from the ambient temperature could cause problems: if the internal temperature is very low, condensate and frosting could be caused by ambient humidity. If the door is often opened under these conditions, the heat exchangers could be obstructed.

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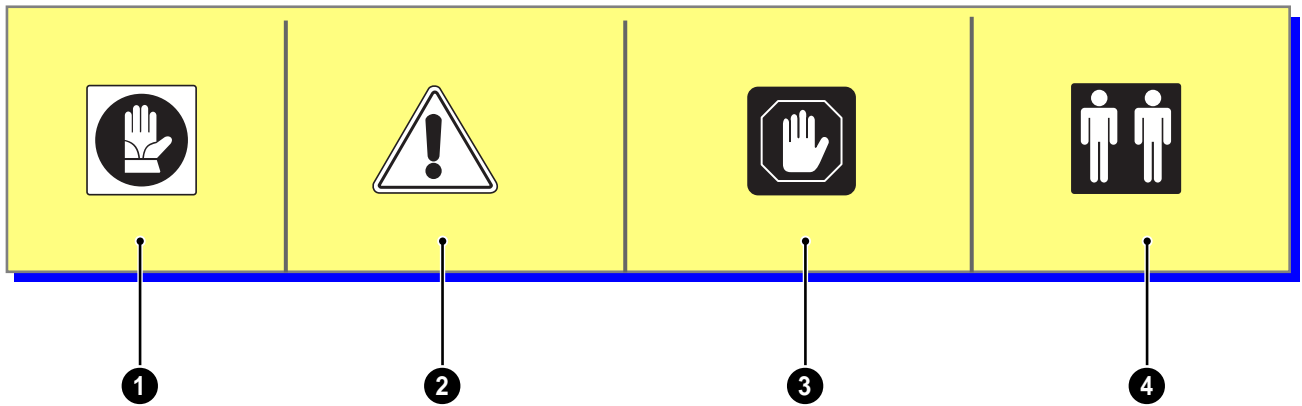## GB 2 WARNINGS

### 2.5 WARNINGS FOR MAINTENANCE

- Always disconnect the machine from the electrical mains before carrying out any maintenance operation.
- Do not use solvents or alcohol to clean the varnished parts as these products could damage the surface.

### 2.6 EXPLANATION OF SYMBOLS

- The symbols shown below may be found on the machine or in this instruction handbook.
- Pay attention to their meaning before going any further.

- 1 **Use gloves**
- 2 **GENERAL DANGER**  
This symbol is always accompanied by an explanation of the danger
- 3 **WARNING**  
Important note
- 4 **Operations that must be carried out by at least two people**

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<sup>(1)</sup> The test is carried out by connecting the thermal couplings (type T) to a mass with the following characteristics:

- material: aluminium alloy
- thermal capacity: 102 J/K
- external surface: 0.009 m<sup>2</sup>


They are positioned at a distance of L/10 from the room walls (L=length of room side), according to IEC 60068-3-5.

Temperature fluctuations over time: maximum difference between the lowest and highest value measured by each sensor (reference value = mean value at the point under study).

Temperature variations inside the space: maximum difference between the mean values measured by the 9 sensors (reference value = set-point value).

<sup>(2)</sup> Measured in front of the source at a distance of 1 m in an environment with no reverberation.

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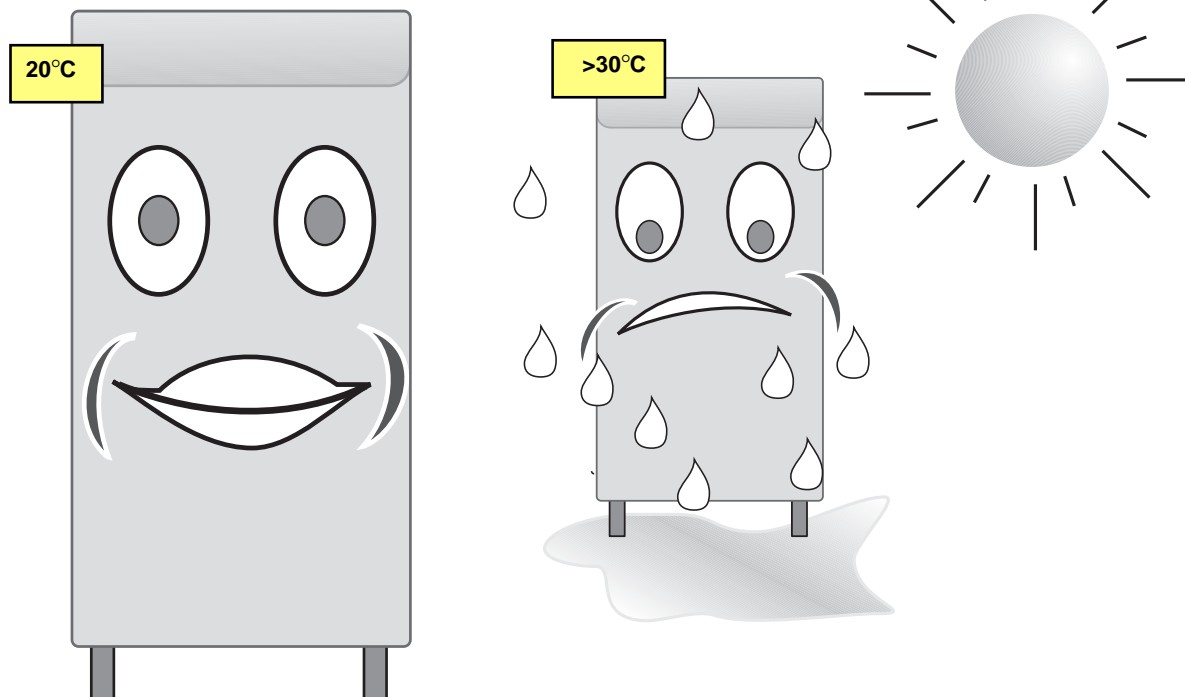#### 3.3 ENVIRONMENTAL CONDITIONS


In order for the appliance to operate correctly, it should be placed in a site with the following requisites:

- far from heat sources,
- far from direct sunlight,
- far from air conditioning systems,
- in a dust-free ambient.

**Ambient temperature** min: +10°C (\*)  
max + 30°C

**Relative humidity RH** max 80%

(\*) The ambient temperature must not be inferior to the chosen set-point.

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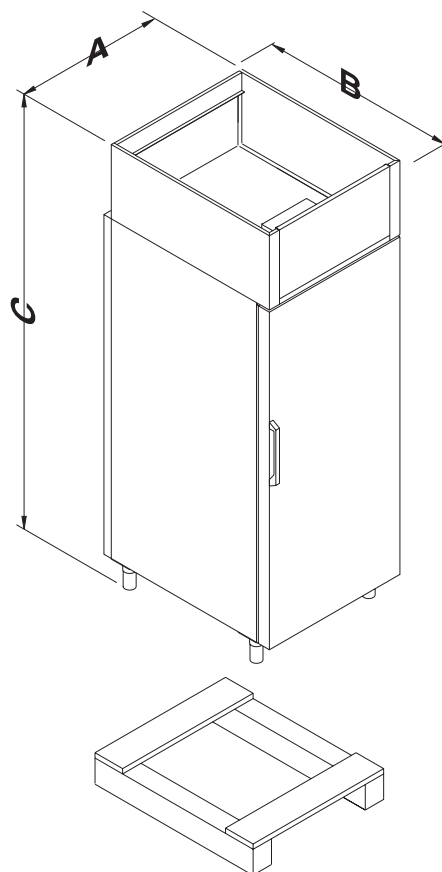## 4 HANDLING AND REMOVAL OF PACKAGING

### 4.1 PERSONNEL REQUISITES

Personnel in charge of handling the machine need no special requisites (just remember the type of packaging).

However, we suggest that this is done by someone who is accustomed to using machines for lifting and transport operations.

### 4.2 MACHINE CONDITIONS

The machine is normally supplied packed and on a pallet. If the machine is delivered by our staff, it may be without packaging. Other types of packaging may be supplied according to the destination and/or customer's needs.

### 4.3 EQUIPMENT NEEDED FOR HANDLING

The machine can be raised and moved by a lift truck, bridge crane, crane or other suitable means with an adequate capacity.

In order to handle the machine, the following rules should be observed:

- Move the machine slowly
- Do not tilt the machine, always keep it in an upright position
- Always make sure you can control the machine during handling

#### Warning:

- Do not turn the machine upside-down
- Do not drag the machine
- Do not shake the machine

#### How to store the machine

It should be kept in a dry environment with an ambient temperature ranging from 0° ÷ 40° C.



Do not place packaged machines one on top of the other and always keep them in an upright position as indicated by the special arrows marked on the packaging itself.

Check that the forks are level and the load is stable before moving the machine.

Make sure that the forks do not hit the machine frame. If it is hit accidentally, check immediately for any damage and if necessary, contact the manufacturer.

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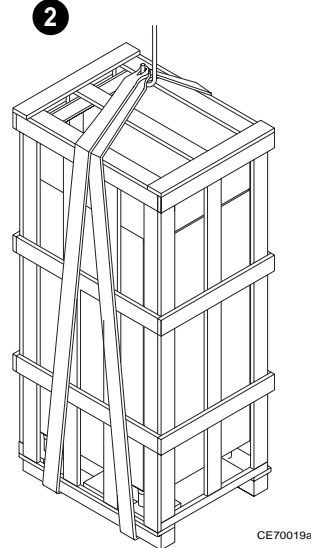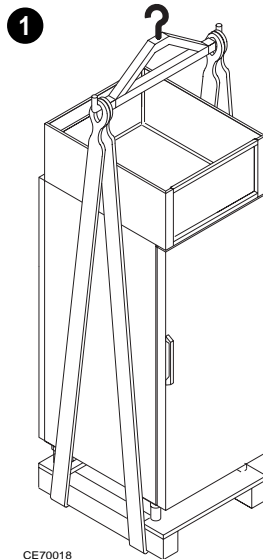## GB 4 HANDLING AND REMOVAL OF PACKAGING

### 4.4 HOW TO RAISE AND MOVE WITH A BRIDGE CRANE

Accessories: textile fibre belts and PVC angle bars.  
Weight: see technical specifications

#### Diagram of slings

- 1 Without packaging
- 2 With packaging

**Do not use place the machine in a sling of metal cables or chains as these could damage it.**

**Never put a sling around a machine without its pallet. In order to prevent the oil in the compressor from flowing into the refrigerating circuit, only transport, stock and handle the equipment in an upright position according to the instructions on the packaging. If the machine is placed on its side, leave it in an upright position for at least 24 hours before switching on.**

### 4.5 REMOVAL OF PACKAGING

The machine is delivered ready-to-go.




Before removing the machine from its packaging, check that it has not been damaged during transport.

If so:

- a) Contact your area representative
- b) Make a written report and if possible attach photographs showing the damage
- c) Send a copy of the report to:
  - Shipping agent insurance company
  - Shipping agent
  - Manufacturer or area representative
- Remove all the packaging and destroy according to the laws in force.
- Eliminate the parts with nails and eliminate plastic bags and plastic sheeting.

**All the packaging materials can be recycled and can be disposed of according to the regulations in force. Take care to eliminate the packaging so that it does not become a danger and throw away the plastic bags, as they could cause children to suffocate.**

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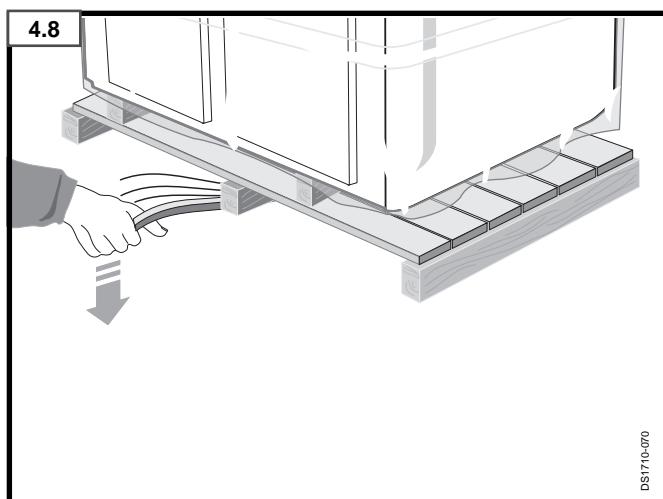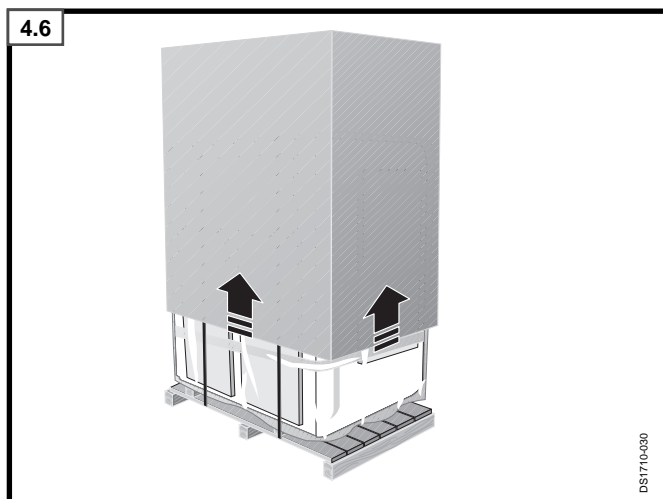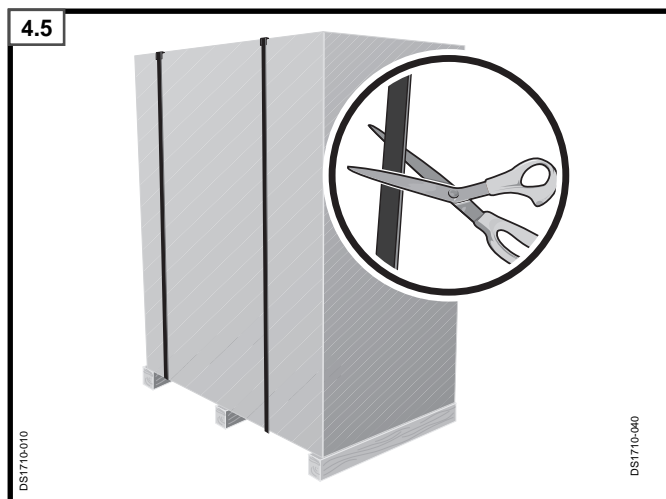## GB 4 HANDLING AND REMOVAL OF PACKAGING

### 4.5.1 How to remove the outer straps (Fig.4.5)

### 4.5.2 How to remove the cardboard casing (Fig.4.6)

### 4.5.3 Removal of plastic materials (Fig.4.7)



When removing the plastic materials, make sure that you don't scratch the surface of the appliance.

Take care not to scratch the surface of the machine.

### 4.5.4 How to remove the wooden side blocks (Fig.4.8)

- Remove the screws that block the ledges on the cabinet.
- Use a special lever to remove the two wooden blocks.

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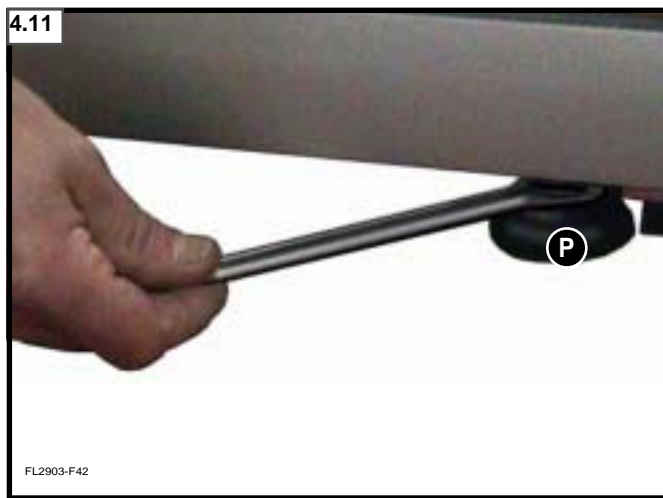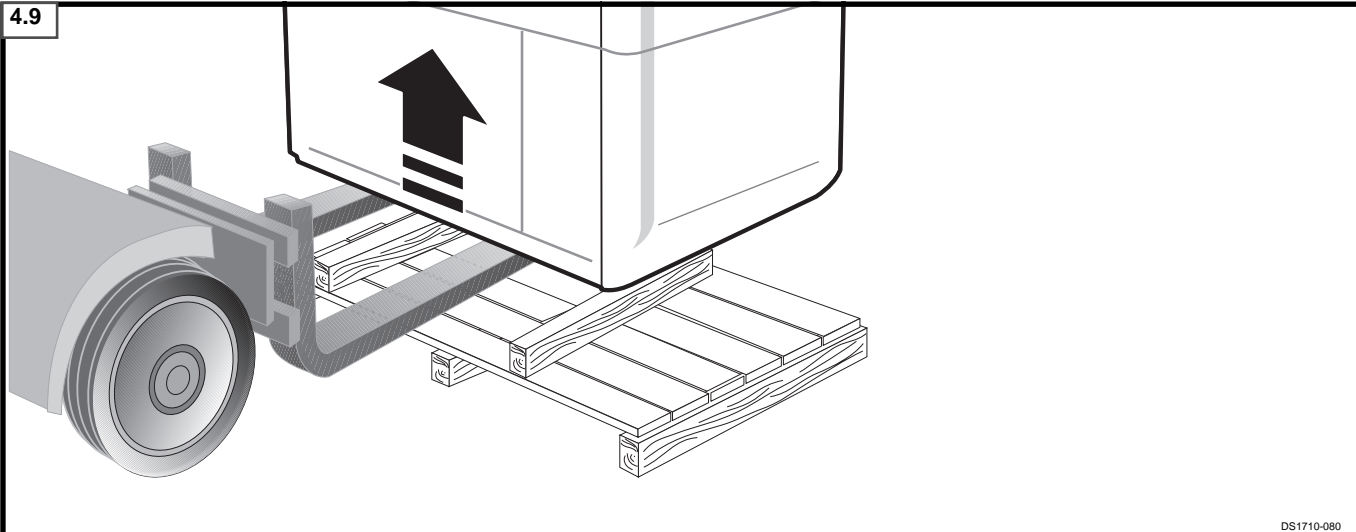## GB 4 HANDLING AND REMOVAL OF PACKAGING

### 4.5.5 How to remove the pallet

- Just raise the machine a few centimetres so that the pallet can be pulled out from underneath.
- The machine can be placed directly on the floor.

**Do not bang or scratch the machine**  
In case of accidental bumps, please check if there is any damage and, should it be necessary, contact the manufacturer.

### MANUAL HANDLING

All models are supported by four spherical wheels and two central rollers; it is therefore possible to move manually.

Before moving the equipment:

- Bear in mind the mass which has to be moved and calculate




the number of people needed to control the mass during movement.

- check that the floor is smooth and completely flat.

### 4.6 HOW TO BLOCK THE MACHINE

The machine must not be fixed to the floor. The machine has been designed to operate on a horizontal floor. In order to prevent the machine from moving accidentally, use the special spanner to adjust the small front feet **P** so that the machine is level and perfectly stable (fig. 4.11).

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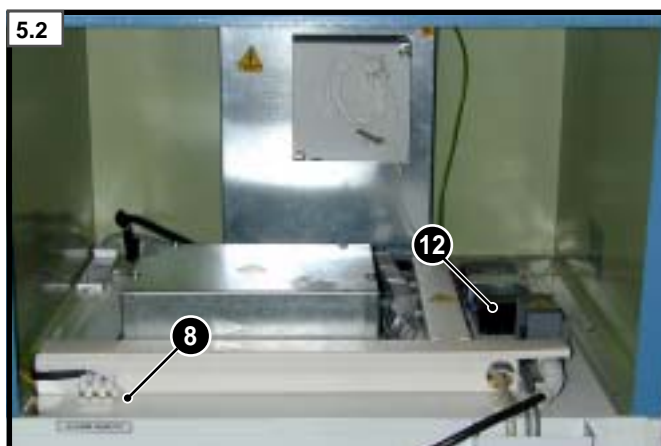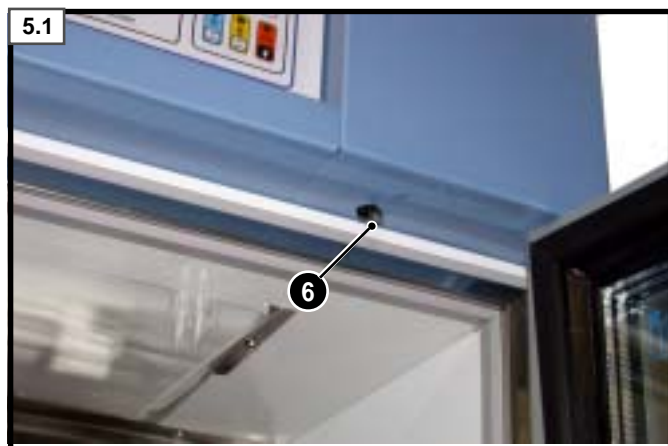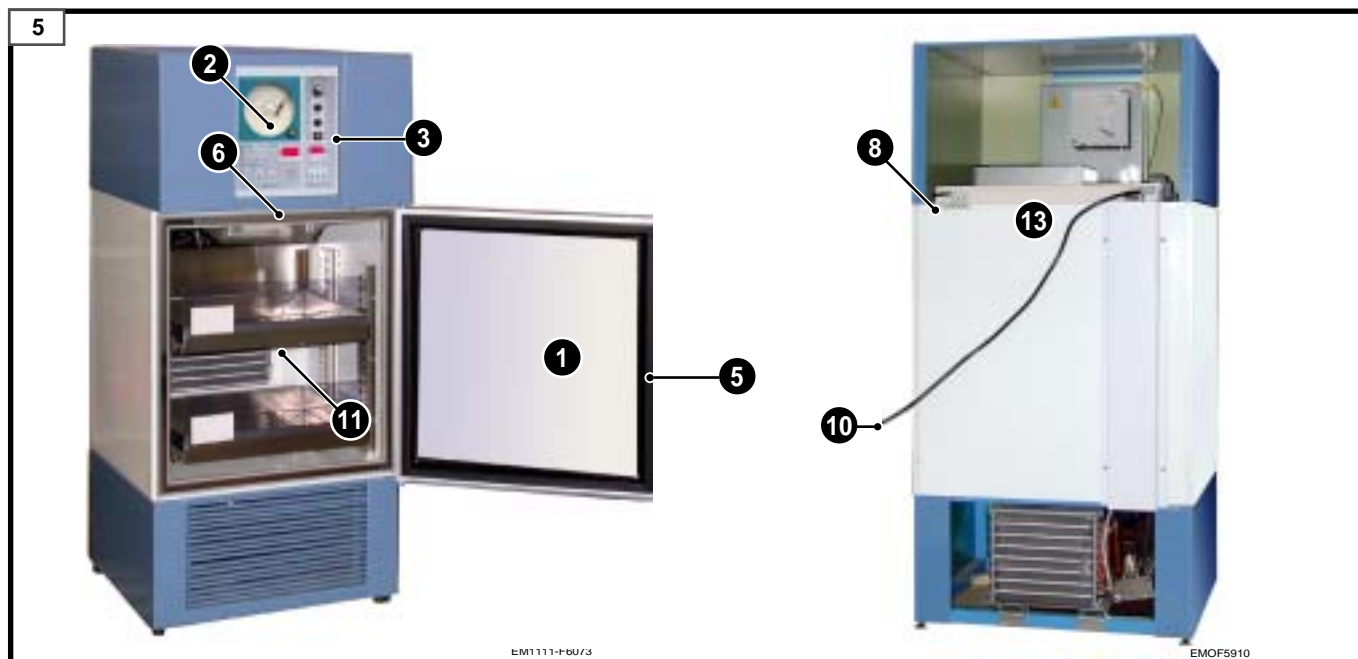## GB 5 DESCRIPTION OF THE SYSTEM

### 5.1 GENERAL VIEW

- 1 Access door/lid
- 2 Chart recorder
- 3 Control panel
- 5 Handle to open door
- 6 Door micro-switch
- 8 Remote alarm connection
- 9 Supply cable connector
- 10 Supply cable
- 11 Drawers
- 12 Buffer battery
- 13 Serial interface RS232 (RS485 optional)

The completely metal structure is made of either pre-painted sheet steel, or, of 304 Aisi stainless steel, on request. The doors are made of anodised steel and triple layered glass (anti-condensate). The structure insulation is obtained with injections of 40kg/m<sup>3</sup> polyurethane foam. The drawers inside the structure can be adjusted in height according to the customer's needs. The cupboard is equipped with internal lighting controlled manually on the front panel.

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## GB 5 DESCRIPTION OF THE SYSTEM

### 5.2 INTERNAL LIGHT

When the door is opened, the door micro-switch **6** automatically switches on the internal light **12**.

#### 12 Internal light

### 5.3 CONTROL SYSTEM

The control panel consists of an electronic board with a micro-processor control and is divided into two sections; one to regulate temperature and defrosting, and the other to run the alarm system and temperature display.

This board is equipped with:

- **5** outputs on relays to run the mechanisms (compressor, defrosting system, evaporator fan, internal refrigerator chamber light, remote control alarm system).
- **3** analog inputs for the PT100 probe (to control temperature, end of defrosting and alarm probe).

A buffer battery (12V 2Ah) is present so as to guarantee the alarm system for approximately 48 hours in the case of a blackout.

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### 6.1 POSITIONING OF THE MACHINE

In order to bring the machine to the installation site, refer to chapter 4.

During installation remember to leave enough space for work and maintenance operations.

Check that the floor is perfectly flat.

Check the ambient conditions of the installation site (see paragr. "3.3 ENVIRONMENTAL CONDITIONS").

**WARNING:**

**The machine has not been designed to operate in areas classed as dangerous (explosive), provided for in the CEI (EIC) norm 64/4. The machine cannot operate in an explosive or corrosive ambient. The equipment being classified as IP 20 is not fit for operation outdoors or in places that are not protected against atmospheric agents.**

Should the user need to place the equipment in places other than those described above, always contact the manufacturer before starting operation.

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## 6.2 ELECTRICAL WIRING

The appliance must be supplied with electricity with the same tension as the one specified on the data tag (see chapter 1 "SUMMARY OF RATING PLATE DATA").

Class 1 insulation, complete with supply cable, without differential switch protection degree  $\leq$  0.5 mA. "THE EQUIPMENT COULD BE A SOURCE OF ELECTRIC DISCHARGE IN THE EVENT OF BREAKDOWN OR INTERRUPTION OF THE PROTECTION CONDUCTOR. THEREFORE, TO IMPROVE SAFETY CONDITIONS, IT COULD BE ALSO INSTALLED (ACCORDING TO THE USER'S REQUIREMENTS) PERMANENTLY" (ref. 6.12.2.1 of CEI/EN 61010-1 adjustments).

This equipment can be considered electrically safe only when it has been correctly wired and an efficient earth system has been installed as foreseen by the safety norms in force.

This fundamental point of safety must be checked and if you are in any doubt, ask for the system to be accurately checked by professional qualified personnel.

In order to ensure that the appliance has been correctly and safely installed, a special, easily accessible socket must be placed near the appliance that is controlled by a multi-polar switch with a 3 mm or more contact opening, in compliance with the safety norms in force, equipped with fuses, with associated differential.


Certain fundamental rules must be followed when using this machine as when using any electrical machine, such as:

- do not touch the equipment with wet or damp hands or feet,
- do not place the machine on normally wet floors,
- do not use extensions in bathrooms or showers,
- do not pull on the supply cable in order to disconnect it from the mains supply.

### 6.2.1 How to connect to the electrical mains supply

**NOTE:** Machine sectioning takes place by disconnecting the plug from the socket.

The machine is supplied with a cable already connected to the panel terminal board; just insert the plug into the mains socket.

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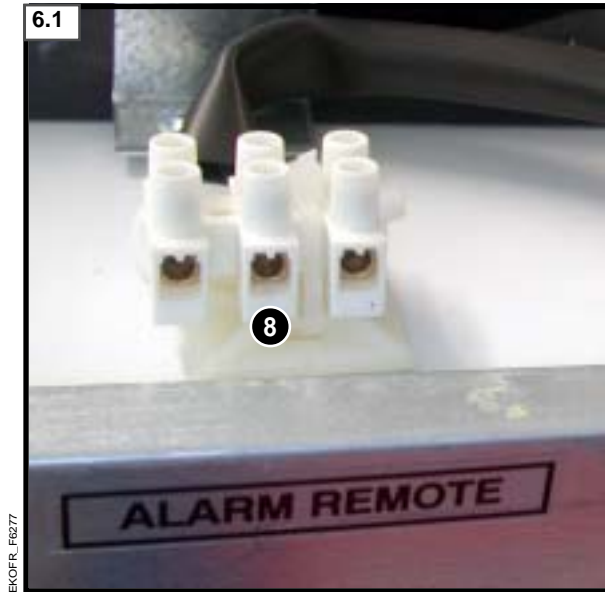## 6 INSTALLATION

### 6.2.2 Supply cable

The supply cable **10** is complete with a **16A** plug in compliance with the adjustments in force (CEI-EN 61010).

The replacement of the supply cable can be made only by qualified personnel of Angelantoni Industrie.

**WARNING! It is strongly recommended to unroll the supply cable for all its length and check that it is not crushed.**

### 6.2.3 Insertion of Remote Alarm

**During this operation the equipment has to be disconnected from the electrical supply. In order to disconnect the equipment from the electrical supply extract the plug from the socket.**

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NB: note the position of terminals **8** at the top of the equipment (Fig.6.1).

### 6.2.4 Serial interface

The serial port (optional) is positioned on the back of the equipment and is indicated by the label **RS232** or **RS485**.

For more information please refer to the attached service handbook of **E2003** control system.

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## GB 7 USE FORESEEN BY THE MANUFACTURER

### 7.1 AIM OF THE MACHINES

These machines are en-bloc refrigerating cupboards, used to conserve diagnostic drugs, biological material, serums and vaccines.

They are furthermore appropriate for the preservation of industrial products in general, within the temperature zones indicated in chapter "3 technical characteristics".

For further information please refer to chapter 3 entitled "Technical specifications".

### 7.2 INTENDED USE OF THE MACHINES

The machines are intended to be used mainly by research laboratories and should only treat the types of materials mentioned in paragraph 7.1.

Before treating any different types of materials, you should consult the manufacturer directly to obtain his specific authorisation.

### 7.3 OPERATOR

No special technical knowledge is needed to use this machine.

### 7.4 USE LIMITS


The machines can be used within the ambient temperature and humidity range described in chapter 3 entitled "TECHNICAL SPECIFICATIONS".

This series of machines must not contain:



- inflammable and/or explosive liquid or solid materials,
- liquid or solid inflammable materials.
- liquid or solid materials which could generate inflammable gases during evaporation or sublimation.
- liquid or solid materials containing chemically aggressive substances,
- liquid or solid materials could generate aggressive fumes during evaporation or sublimation.

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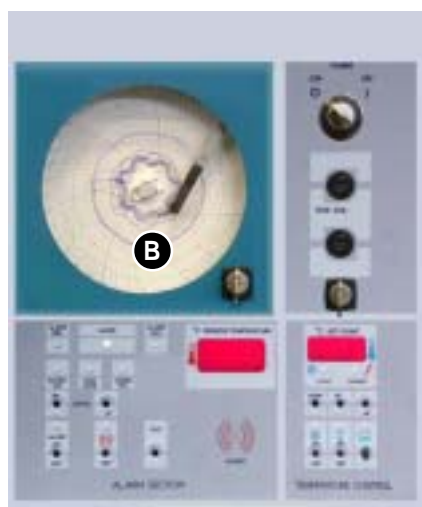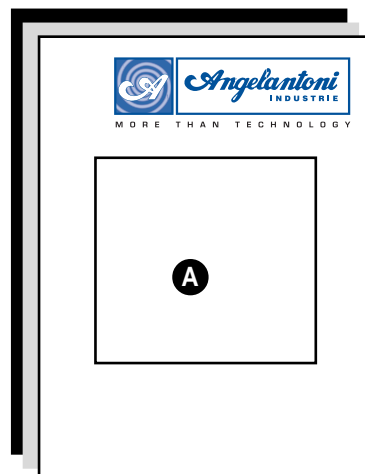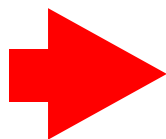
## GB 8 START-UP

- Check that the mains supply voltage is the same as the one marked on the rating plate.

### 8.1 E2003 SYSTEM

*In order to set and start up the chamber with the E2003 system, please refer to the attached instruction handbook.*

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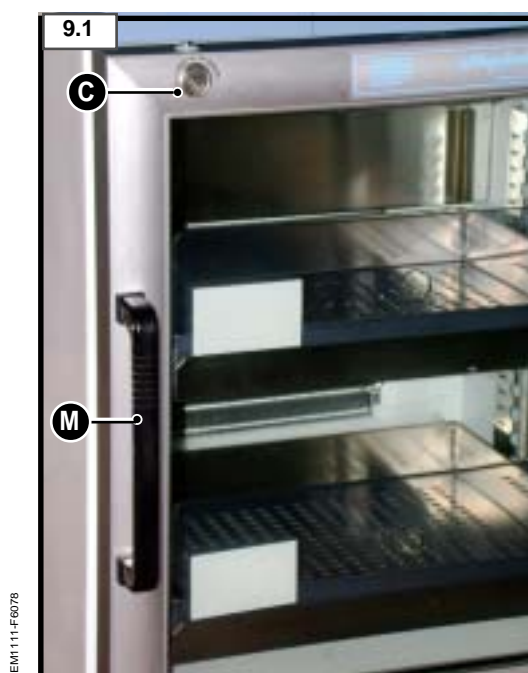## GB 9 USE

- In order to gain access inside the chamber, turn key **C** and pull handle **M** (fig. 9.2).

### 9.1 PRODUCT LOADING

Take care to fasten the shelves carefully and do not place the containers directly on the floor of the refrigerating chamber. Insert the objects into the refrigerating chamber a few at a time once the temperature has become stabilised.

Always leave spaces between the stored materials and do not obstruct the air flow.

Ambient temperature material (e.g. +25°C) should be put inside the freezer gradually and in small quantities to ensure proper storage; allow the equipment to freeze the products inside the refrigerating compartment before introducing new material to be cooled.

If this important rule is not observed, the temperature inside the freezer could suddenly increase and cause previously stored products to deteriorate.

When material is put in or taken out, it must be done very quickly because each time the door is opened, there will be continual unintentional changes in the air so that the temperature will rise and a large amount of ice will form on the cooler.


**Do not store products imbued with liquid substances inside the refrigerator since they could produce vapours during the first cooling phase, leading to ice formation on evaporator.**

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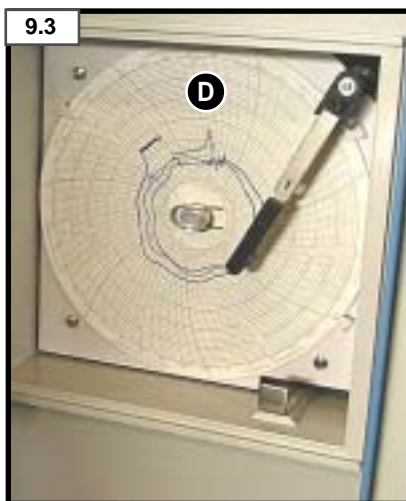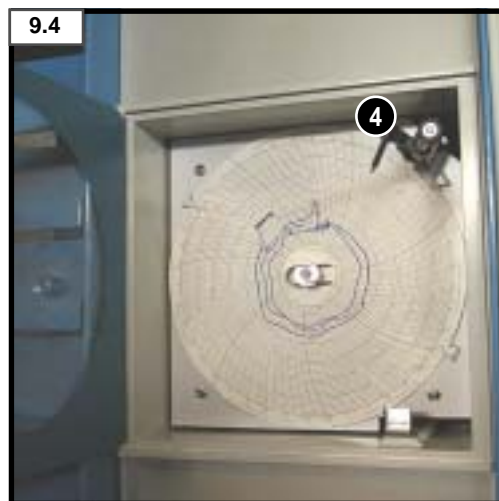
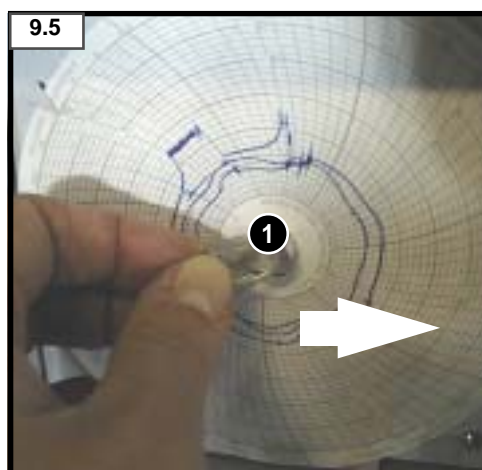
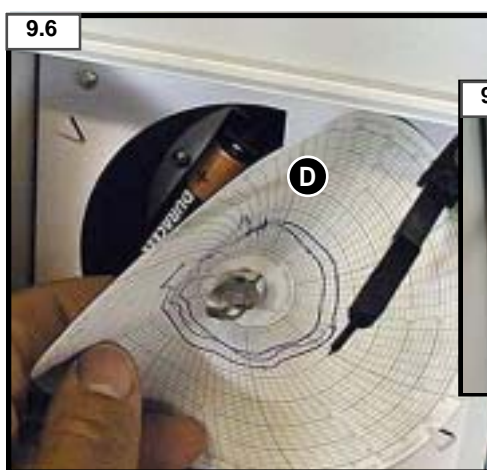
## GB 9 USE

### 9.2 TEMPERATURE CHART RECORDER

Freezers are supplied with a temperature chart recorder **2** upon customer request.

The standard 7-day battery-operated recorder comes complete with a clock run by a 1.5V high-capacity accumulator lasting approximately 6 months.

The temperature chart recorder is moreover supplied complete with no. 1 printing cartridge and no. 55 chart disks ensuring 1-year coverage.

Contact Angelantoni's Assistance Service for additional printing cartridges and chart disks.

#### 9.2.1 Use of the recorder

The recorder has a battery powered clock.

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
- Lift the pen **4**.
- Move the lock **1** in the direction shown by the arrow and lift it.
- Remove the chart disk and install the new chart disk in the central pin. Move the lock **1** to the previous position and lower the pen **4**.

#### Battery substitution (every 6 months)

- Follow the steps described in points **9.4**, **9.5**, **9.6** and move the diagram disk without removing it.
- Replace the battery when necessary with a new 1.5V (Mod. AA LR6) battery.

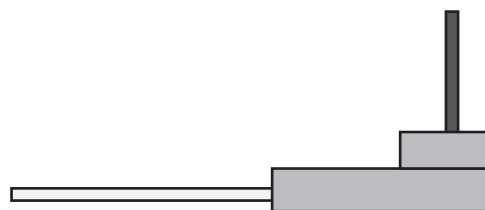**PRIMA PENNA**  
**FIRST PEN**
**SECONDA PENNA**  
**SECOND PEN**
**MEDIO**  
**MEDIUM**
**LUNGO**  
**LONG**

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**GB 9 USE**
**WRITING**
**Writing cartridges.**

The instruments are supplied with pen points already inserted into their respective rods. In order to start up, proceed as follows:

- remove the protective, plastic caps from the cartridge writing tips.

**Writing cartridges replacement**

If the writing should appear faulty, replace the pen points.

Proceed as follows:





- Hold the rod firmly and remove the pen point that has to be replaced by pulling it downwards.

- Insert a new pen point of the same colour and the same length.

**To order the nib, please specify the code as indicated in the following table.**

**The code number to order the recorder disk is printed on the centre of the recorder disk itself.**

RECORDER TYPE	COLOURS	
	FIRST PEN	SECOND PEN
<b>SINGLE PEN</b>	MEDIUM BLUE code N° 133962	
<b>TWO PENS</b>	MEDIUM RED code N° 133969	LONG BLUE code N° 133968

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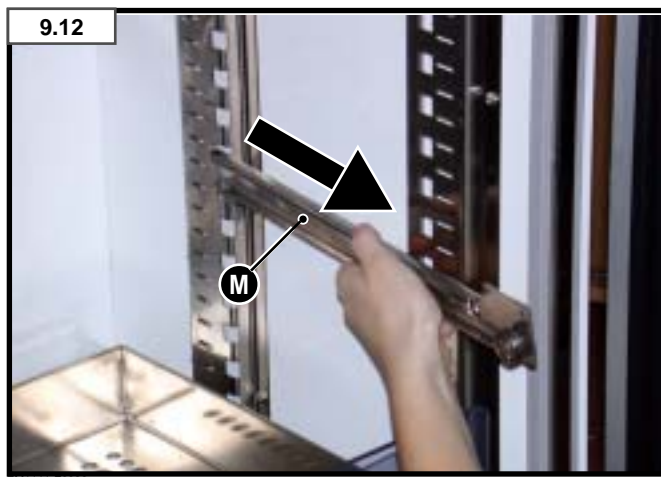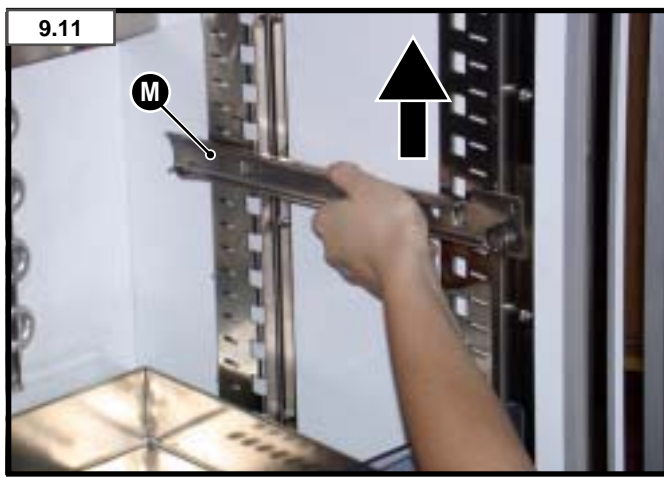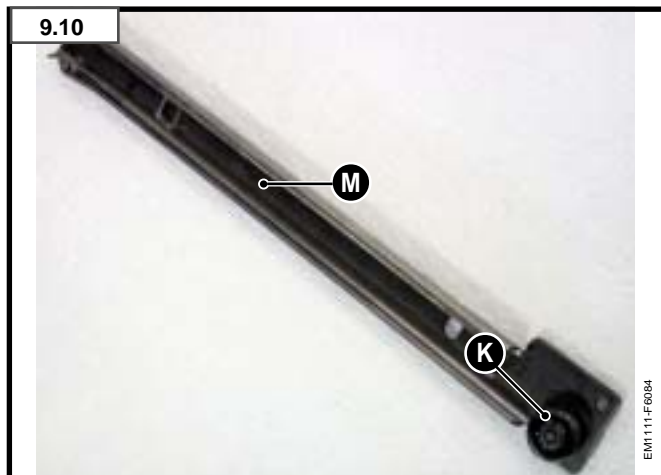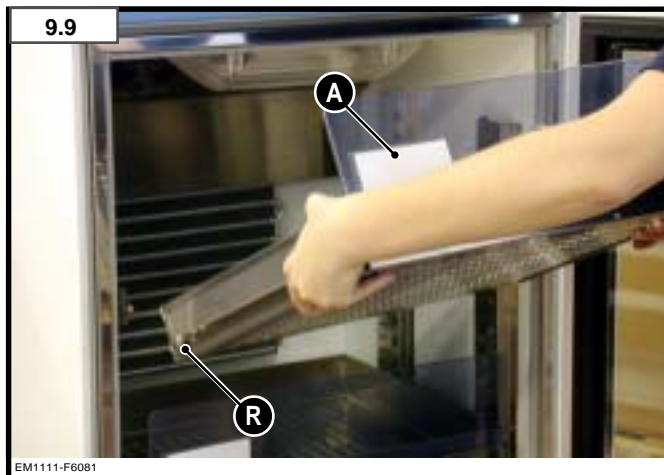## GB 9 USE

### 9.3 HOW TO INSERT AND CHANGE THE DRAWERS

The drawer **A** runner is fixed. The wheel **R** is installed at the end of the runner.

The runners **M** can be moved in vertical position.

- Extract the drawers (Fig. 9.9).

#### Extraction of the runners (Fig. 9.11):

- Lift the frontal part of the runner.
- Extract the back part of the runner (Fig.9.12).
- Re-mount the runners carrying out these operations in reverse.
- Re-mount the drawers.



### • CHANGE THE DRAWERS

**WARNING!** If the machine has been loaded (i.e. there is material in the drawers) open one drawer at a time. Do not open all the drawers at the same time, otherwise the excess weight may unbalance the machine and cause it to fall over face down.

### 9.4 CONNECTION OF OPTIONAL PT100 PROBE (optional)

When this optional is provided, on the back of the equipment there is an electrical box with the indication "PT100 PROBE". Inside this box there are the terminals of the PT100 wires (the sensing part of the probe is positioned inside the equipment). In order to connect the probe to Your system, open the box with a screwdriver and make the connection.

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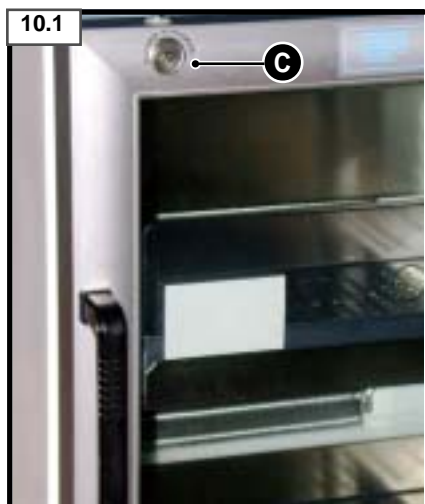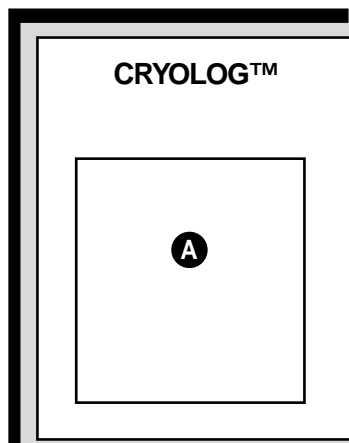## GB 9 USE

### 9.5 DATA DOWNLOADING VIA SERIAL PORT (optional)

- Connect your computer to the serial port **13** (see paragr. 6.2.4).
- Switch on your computer and start the programme.
- In order to set and start up the software, please refer to the attached instruction handbook according to the configuration of your machine.


**A :** Cryolog™ software instruction handbook

## GB 10 SAFETY DEVICES- CHECK AND SET-UP

SAFETY DEVICE	CHECK	SET-UP
Key <b>C</b> to lock handle (Fig.10.1)	This key ensures that the machine cannot be opened either accidentally or on purpose by unauthorized personnel. The machine may be at very low temperatures.	<ul style="list-style-type: none"> <li>• Close the door, remove the key and check that door will not open.</li> </ul>
Fuse thermostat	The fuse thermostat is set at the factory at a temperature above the chamber work temperature. When it triggers, it stops the machine. In order to reset it, it has to be replaced by the technical assistance service.	<ul style="list-style-type: none"> <li>• If the machine switches on, this shows that the fuse thermostat is working correctly.</li> </ul>
Door micro-switch (Fig.10.2)	When the door is opened, the micro-switch <b>6</b> automatically turns on the internal light and turns off evaporator ventilation.	<ul style="list-style-type: none"> <li>• Close the door.</li> </ul>

**Other safety systems concerning the cooling system cannot be checked or set by the user and are, therefore, described in chapters 11 and 12.**

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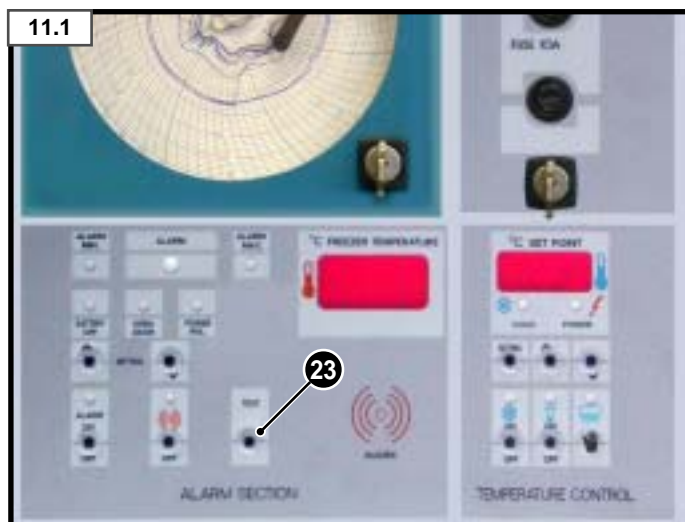## GB 11 MAINTENANCE

Maintenance must be carried out by qualified personnel.

Some maintenance operations may be carried out when the machine is working and therefore all the necessary safety precautions should be taken.

With the exception of specific cases maintenance operations must be carried out in the following conditions:

- machine switched off
- chamber temperature the same as ambient temperature
- machine disconnected from the mains supply
- supply valves or taps above the machine closed.

### 11.1 ALARM SYSTEM CONTROL

Please refer to the attached instruction handbook E2003 system. (ref. 23 - TEST ALARM).


### 11.2 DEFROSTING

Please refer to the attached instruction handbook E2003 system.

### 11.3 HOW TO CLEAN THE STRUCTURE

During this operation the equipment has to be disconnected from the electrical supply. In order to disconnect the equipment from the electrical supply extract the plug from the socket.

- Remove the shelves.
- Clean with a damp cloth dipped in a neutral detergent and then use another damp cloth dipped in pure water.
- Replace the shelves.

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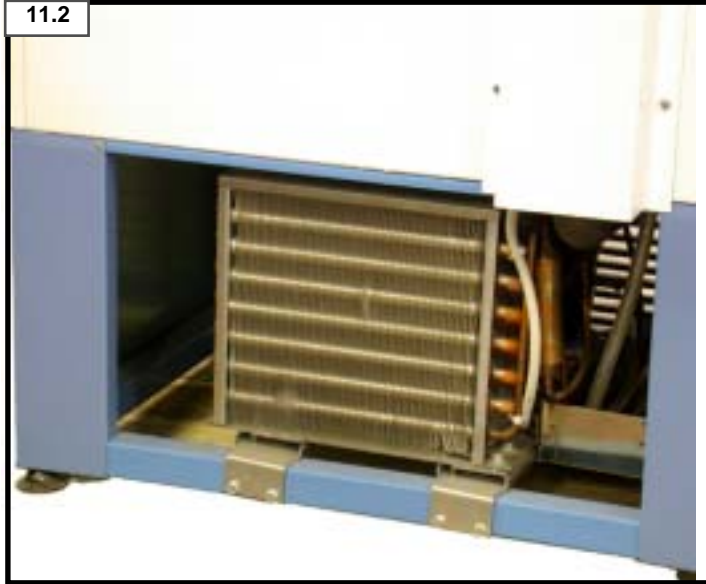## 11 MAINTENANCE

## 11.3.1 How to clean the gaskets

During this operation the equipment has to be disconnected from the electrical supply. In order to disconnect the equipment from the electrical supply extract the plug from the socket.

If the gaskets are dirty, they will not function correctly.

Clean them with a damp cloth dipped in a neutral detergent and then use another damp cloth dipped in pure water.

Never spray any part of the equipment with water as this could damage the electrical insulation and prevent the machine from operating correctly.

## 11.3.2 How to clean the condenser

During this operation the equipment has to be disconnected from the electrical supply. In order to disconnect the equipment from the electrical supply extract the plug from the socket.

Use a soft brush or jets of compressed air to clean the condenser once every 6 months. Models with two compartments have two condensers (one for each compartment).

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## 11.4 COOLING SYSTEM MAINTENANCE

Each machine is equipped with a refrigerating system using a hermetically sealed compressor, and therefore, (if there are any accidental leaks) no refrigerant or non-freezing oil has to be added.

We advise you to have it checked once a year by qualified personnel.

## 11.5 MAINTENANCE OF THE ELECTRICAL-ELECTRONIC PART

The electronic control is a very simple, reliable, integrated system in a single container, so no maintenance is required.

We advise you to have it checked once a year by qualified personnel in order to see that the alarm system is operating correctly, that the visualised temperature remains linear, and to replace the compressor control relay (every 5 years).

Our company provides its own economic, subscriber "SPECIALIZED TECHNICAL ASSISTANCE SERVICE".

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
The maintenance operations authorizes are those indicated by the producer are indicated in the “Regular maintenance” chapter. The information given in the column “SOLUTION” in the tables below does not authorize any operations if these compromise safety; the information is given to help any specialized technicians find the fault.

PROBLEM	PROBABLE REASON	SOLUTION
The machine does not start.	<ul style="list-style-type: none"> <li>There is no supply from the mains.</li> <li>The main switch is in “OFF” position</li> <li>The plug is not connected to the main</li> <li>Some other reason not mentioned above</li> </ul>	<ul style="list-style-type: none"> <li>Check and proceed accordingly</li> <li>Turn ON the main switch</li> <li>Connect the plug to the main</li> <li>Call your technical assistance service.</li> </ul>
The equipment does not reach the set temperature	<ul style="list-style-type: none"> <li>Setpoint programming error.</li> <li>There is too much material in the test chamber.</li> <li>The test compartment contain too much and/or too wet material.</li> <li>Malfunctioning of the regulation system.</li> <li>Malfunctioning of the cooling system.</li> </ul>	<ul style="list-style-type: none"> <li>Check the set temperature value following the procedure as in (ref. par.8.1).</li> <li>Reduce the quantity.</li> <li>Allow temperature values to stabilise.</li> <li>Call the technical assistance service.</li> <li>Call the technical assistance service.</li> </ul>
The equipment is too noisy.	<ul style="list-style-type: none"> <li>The equipment is not evenly laid on the floor.</li> <li>Contact with foreign bodies.</li> <li>Damaged mechanical parts.</li> </ul>	<ul style="list-style-type: none"> <li>Remove the cause.</li> <li>Check that foreign bodies such as boxes or other structures do not come into contact with the equipment so as to avoid vibrations during functioning.</li> <li>Check that moving mechanical parts such as fans are not sources of noise.</li> </ul>
The display shows alarm codes.	<ul style="list-style-type: none"> <li>Refer to paragraph (ref. 8.1) of these operating instructions.</li> </ul>	<ul style="list-style-type: none"> <li>Call the technical assistance service and specify the type of code displayed as well as the possible cause.</li> </ul>
The light does not switch on.	<ul style="list-style-type: none"> <li>The lamp is worn out or burnt out.</li> </ul>	<ul style="list-style-type: none"> <li>Replace the lamp.</li> </ul>

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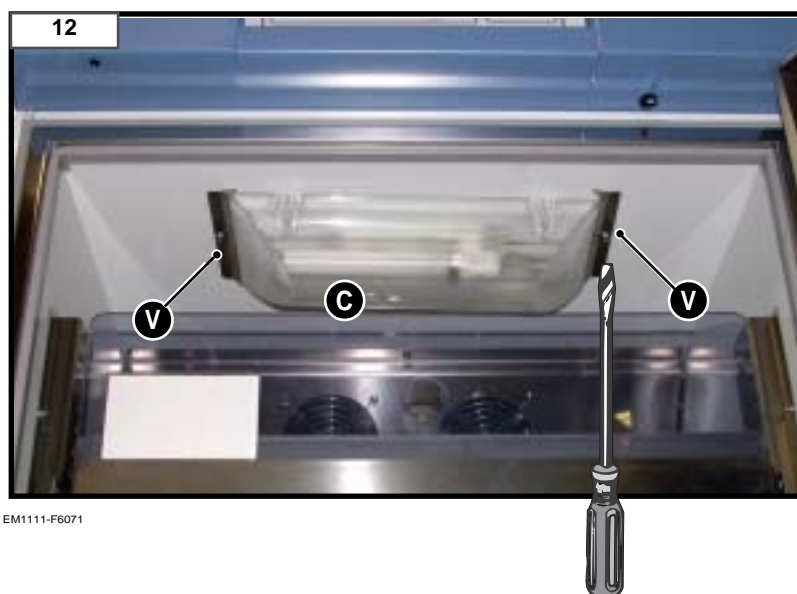## GB 12 TROUBLESHOOTING

### 12.1 HOW TO REPLACE THE LAMP

During this operation the equipment has to be disconnected from the electrical supply. In order to disconnect the equipment from the electrical supply extract the plug from the socket.

**WARNING!** If the machine has been loaded (i.e. there is material in the drawers) open one drawer at a time. Do not open all the drawers at the same time, otherwise the excess weight may unbalance the machine and cause it to fall over face down.

- Remove screws **V** and then remove cover **C**.
- Remove the lamp from its holder and replace it with a new one.
- Replace cover **C** and tighten screws **V**.

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The machine can be disassembled in order to move it from one site to another or in order to be scrapped. Transport must be carried out in compliance with the norms described in this handbook in the chapters "HANDLING" and "INSTALLATION".

### 13.1 DISASSEMBLY

- In order to disconnect the machine from the mains supply, follow the procedure described for the electrical wiring in reverse.
- Clean the surfaces accurately, oil the parts that could rust and cover the machine so as to protect it from dust, dirt and humidity.

### 13.2 SCRAPPING

- In order to disconnect the machine from the mains supply, follow the procedure described for the electrical wiring in reverse.
- Discharge all the refrigerant from the refrigerating circuit.
- Remove all the oil from the motors.


**WARNING! This product should not be left in the surrounding environment:**

- refrigerating substances;
- the compressor oil;
- polyurethane.

**It would be advisable to contact specialised companies for the disposal and the recycling of the above materials.**

- Divide the remaining parts according to their type and destroy in compliance with the laws in force.
- To remove and transport the machine please refer to the procedure described in the chapters "HANDLING" and "INSTALLATION".

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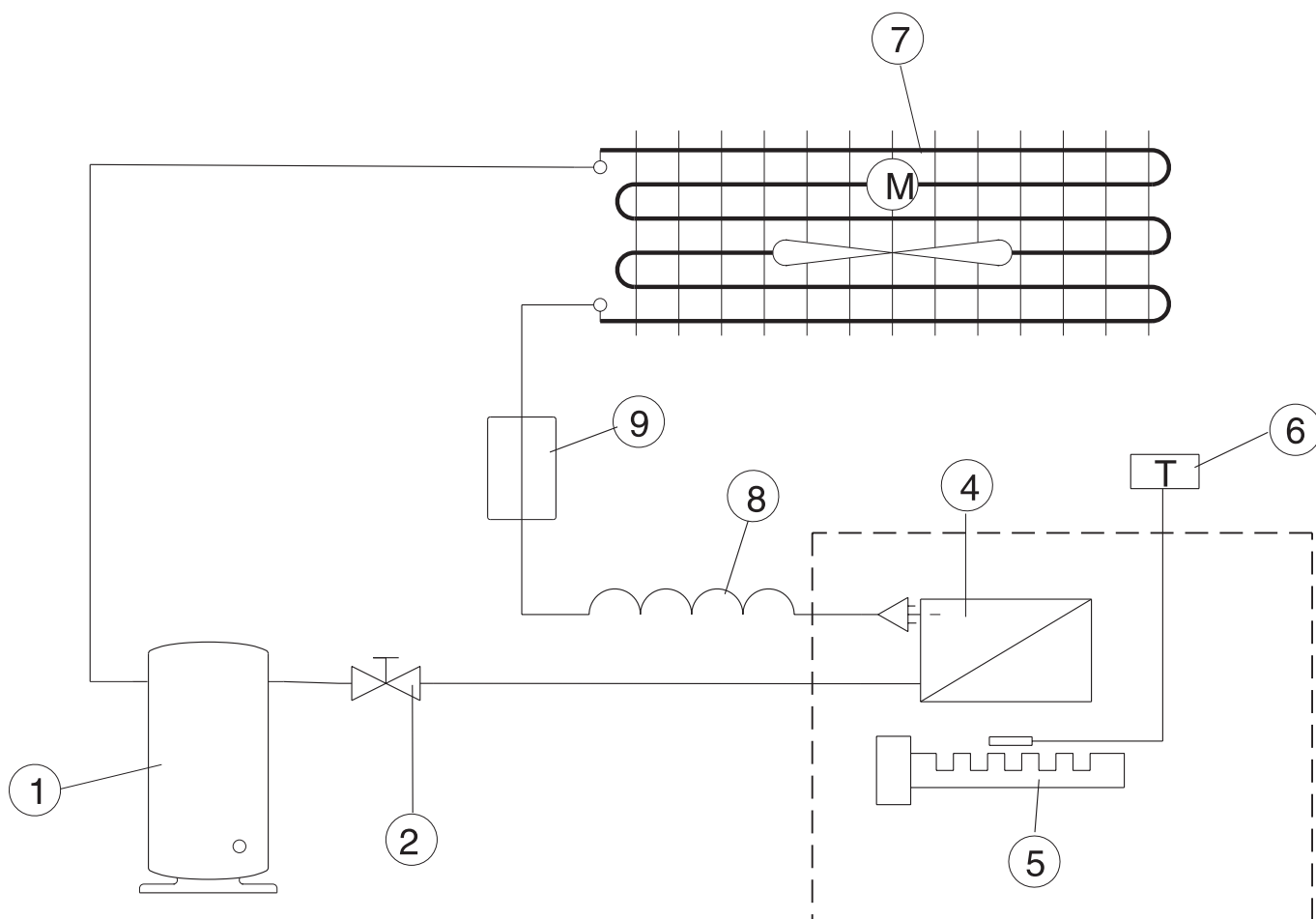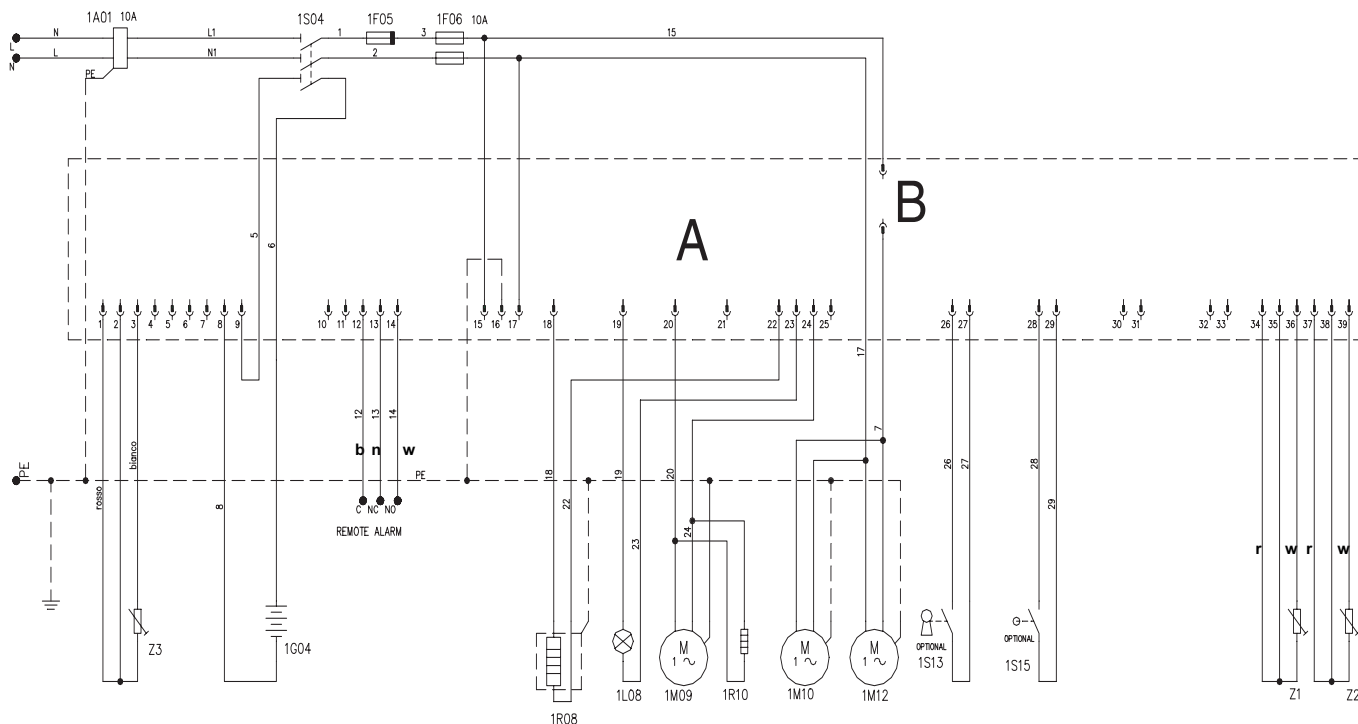
## GB 14 DIAGRAMS

### 14.2 ELECTRIC DIAGRAM

<b>A</b>	Board E2003	<b>1R10</b>	Door frames heaters
<b>B</b>	Compressor relay	<b>1S04</b>	Main switch
<b>1A01</b>	Antijamming filter	<b>1S13</b>	Control block switch
<b>1F05</b>	Thermal fuse	<b>1S15</b>	Door micro switch
<b>1F06</b>	Line fuses	<b>Z1</b>	Regulation probe
<b>1G04</b>	Buffer battery	<b>Z2</b>	Evaporator probe
<b>1L08</b>	Internal light	<b>Z3</b>	Alarm probe
<b>1M09</b>	Evaporator ventilation motor	<b>b</b>	Blue cable
<b>1M10</b>	Condenser ventilation motor	<b>n</b>	Black cable
<b>1M12</b>	Compressor motor	<b>w</b>	White cable
<b>1R08</b>	Defrosting heaters	<b>r</b>	Red cable

## DET NORSKE VERITAS

### QUALITY MANAGEMENT SYSTEM CERTIFICATE

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This certificate is valid for the following products or services:

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Design, manufacture, sales and after sales services of simulated environmental test chambers and climatic plants, refrigerating equipment for biomedical applications and scientific research, cooling plants for industrial processes

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Lead Auditor: **Flavio Franco**

Settore EA: 18

**SINCERT**  
ACCREDITED ORGANISM CERTIFICATION

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**Leonardo Omodeo Zorini**  
Management Representative

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## DET NORSKE VERITAS

### ENVIRONMENTAL MANAGEMENT SYSTEM CERTIFICATE

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CONFORMS TO THE ENVIRONMENTAL MANAGEMENT SYSTEM STANDARD

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Certification has been granted in conformity with the SINCERT Technical Regulation RT - 09

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Design, manufacture, sales and after sales services of simulated environmental test chamber and climatic plants, refrigerating equipment for biomedical applications and scientific research, cooling plants for industrial processes by cutting, punching, plying, welding, mechanical assembling, electrical assembling, refrigerators equipment assembling and testing

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**Vitore Marangon**  
Management Representative

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**Det Norske Veritas Italia S.r.l.**